

Zyto Dot 2C SPEC USP6 Break Apart Probe

REF C-3077-100



10 (0.1 ml)

For the qualitative detection of translocations involving the human USP6 gene at 17p13.2 by chromogenic *in situ* hybridization (CISH)



In vitro diagnostic medical device according to EU directive 98/79/EC

1. Intended use

The <u>Zyto Dot 2C SPEC USP6 Break Apart Probe</u> (**PD56**) is intended to be used for the qualitative detection of translocations involving the human USP6 gene at 17p13.2 in formalin-fixed, paraffin-embedded specimens by chromogenic *in situ* hybridization (CISH). The probe is intended to be used in combination with the <u>Zyto Dot 2C CISH Implementation Kit</u> (Prod. No. C-3044-10/-40).

Interpretation of the results must be made within the context of the patient's clinical history with respect to further clinical and pathologic data of the patient by a qualified pathologist.

2. Clinical relevance

Translocations affecting USP6 (ubiquitin specific peptidase 6) have been initially found in primary aneurysmal bone cysts (ABC), a benign, but locally aggressive bone lesion that occurs predominantly during the first two decades of life. USP6 rearrangements are restricted to spindle cells in primary ABC, indistinguishable from surrounding normal spindle cells. The resulting fusion genes detected are formed by juxtaposition of the USP6 coding sequences to the highly active promoter sequences of several partner genes, as e.g. CDH11, COL1A1, OMD, TRAP150, and ZNF9, leading to the transcriptional upregulation of USP6. No true fusion genes are formed. More recently, nodular fasciitis (NF), another mesenchymal lesion, has been tested positive for USP6 rearrangements. NF is a subcutaneous pseudosarcomatous myofibroblastic proliferation of unknown pathogenesis that regresses spontaneously when not surgically resected. The translocation results in the fusion of the promoter region of MYH9 located on 22q12.3 to the entire coding sequence of USP6 and subsequently in upregulated USP6 expression. For both lesions it is assumed that the detection of USP6 rearrangements by chromogenic in situ hybridization might represent a valuable diagnostic tool.

3. Test principle

The chromogenic *in situ* hybridization (CISH) technique allows the detection and visualization of specific nucleic acid sequences in cell preparations. Hapten-labeled nucleotide fragments, so called CISH probes, and their complementary target sequences in the preparations are co-denatured and subsequently allowed to anneal during hybridization.

Afterwards, unspecific and unbound probe fragments are removed by stringency washing steps. Duplex formation of the labeled probe can be visualized using primary (unmarked) antibodies, which are detected by secondary polymerized enzyme-conjugated antibodies. The enzymatic reaction with chromogenic substrates leads to the formation of colored precipitates. After counterstaining the nucleus with a nuclear dye, hybridized probe fragments are visualized by light microscopy.

4. Reagents provided

The Zyto Dot 2C SPEC USP6 Break Apart Probe is composed of:

- Digoxigenin-labeled polynucleotides (~0.50 ng/µl), which target sequences mapping in 17p13.2* (chr17:4,825,753-5,017,582) distal to the USP6 breakpoint region (see Fig. 1).
- Dinitrophenyl-labeled polynucleotides (~0.75 ng/μl), which target sequences mapping in 17p13.2* (chr17:5,087,046-5,361,104) proximal to the USP6 breakpoint region (see Fig. 1).
- Formamide based hybridization buffer

^{*}according to Human Genome Assembly GRCh37/hg19

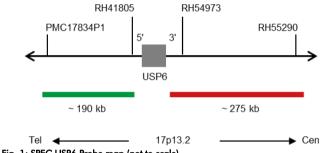


Fig. 1: SPEC USP6 Probe map (not to scale)

The Zyto Dot 2C SPEC USP6 Break Apart Probe is available in one size:

C-3077-100: 0.1 ml (10 reactions of 10 μl each)

5. Materials required but not provided

- Zyto Dot 2C CISH Implementation Kit (Prod. No. C-3044-10/-40)
- Positive and negative control specimens
- Microscope slides, positively charged
- Water bath (80°C, 98°C)
- Hybridizer or hot plate
- Hybridizer or humidity chamber in hybridization oven
- Adjustable pipettes (10 μ l, 1000 μ l)
- Staining jars or baths
- Timer
- Calibrated thermometer
- Ethanol or reagent alcohol
- Xylene
- Methanol 100%
- Hydrogen peroxide (H₂O₂) 30%
- Deionized or distilled water
- Coverslips (22 mm x 22 mm, 24 mm x 32 mm)
- Rubber cement, e.g., <u>Fixogum Rubber Cement</u> (Prod. No. E-4005-50/-125) or similar
- Adequately maintained light microscope (400-630x)

6. Storage and handling

Store at 2-8°C in an upright position. Return to storage conditions immediately after use. Do not use reagents beyond expiry date indicated on the label. The product is stable until expiry date indicated on the label when handled accordingly.

7. Warnings and precautions

- Read the instructions for use prior to use!
- Do not use the reagents after the expiry date has been reached!
- This product contains substances (in low concentrations and volumes) that are harmful to health and potentially infectious. Avoid any direct contact with the reagents. Take appropriate protective measures (use disposable gloves, protective glasses, and lab garments)!

Vers. 3.1.1 EN

- Report any serious incident that has occurred in relation to the product to the manufacturer and the competent authority according to local regulations!
- If reagents come into contact with skin, rinse skin immediately with copious amounts of water!
- A material safety data sheet is available on our homepage (www.zytovision.com).
- Do not reuse reagents, unless reuse is explicitly permitted!
- Avoid any cross-contamination and micro-bacterial contamination of the reagents!
- The specimens must not be allowed to dry during the hybridization and washing steps!

Hazard and precautionary statements:

The hazard-determining component is formamide.



Danger

H351	Suspected of causing cancer.
H360FD	May damage fertility. May damage the unborn child.
H373	May cause damage to organs through prolonged or repeated exposure.
P201	Obtain special instructions before use.
P202	Do not handle until all safety precautions have been read and understood.
P260	Do not breathe dust/fume/gas/mist/vapours/spray.
P280	Wear protective gloves/protective clothing/eye protection/face protection.
P308+P313	IF exposed or concerned: Get medical advice/attention.
P405	Store locked up.

8. Limitations

- For in vitro diagnostic use.
- For professional use only.
- For non-automated use only.
- The clinical interpretation of any positive staining, or its absence, must be done within the context of clinical history, morphology, other histopathological criteria as well as other diagnostic tests. It is the responsibility of a qualified pathologist to be familiar with the CISH probes, reagents, diagnostic panels, and methods used to produce the stained preparation. Staining must be performed in a certified, licensed laboratory under the supervision of a pathologist who is responsible for reviewing the stained slides and assuring the adequacy of positive and negative controls.
- Specimen staining, especially signal intensity and background staining, is dependent on the handling and processing of the specimen prior to staining. Improper fixation, freezing, thawing, washing, drying, heating, sectioning, or contamination with other specimens or fluids may produce artefacts or false results. Inconsistent results may result from variations in fixation and embedding methods, as well as from inherent irregularities within the specimen.
- The probe should be used only for detecting loci described in chapter 4. "Reagents provided".
- The performance was validated using the procedures described in these instructions for use. Modifications to these procedures might alter the performance and have to be validated by the user.

9. Interfering substances

The following fixatives are incompatible with ISH:

- Bouin's fixative
- B5 fixative
- Acidic fixatives (e.g., picric acid)
- Zenker's fixative
- Alcohols (when used alone)
- Mercuric chloride
- Formaldehyde/zinc fixative
- Hollande's fixative
- Non-buffered formalin

10. Preparation of specimens

Recommendations:

- Avoid cross-contamination of samples in any step of preparation as this may lead to erroneous results.
- Fixation in 10% neutrally buffered formalin for 24 h at room temperature (18-25°C).
- Sample size ≤ 0.5 cm³.
- Use premium quality paraffin.
- Embedding should be carried out at temperatures lower than 65°C.
- Prepare 3-5 μ m microtome sections.
- Use positively charged microscope slides.
- Fix tissue sections for 2-16 h at 50-60°C.

11. Preparatory treatment of the device

The product is ready-to-use. No reconstitution, mixing, or dilution is required. Bring probe to room temperature (18-25°C) and mix briefly before use.

12. Assay procedure

Specimen pretreatment

Perform specimen pretreatment (e.g., dewaxing, proteolysis) according to the instructions for use of the <u>ZytoDot 2C CISH Implementation Kit</u>.

Denaturation and hybridization

- 1. Pipette $10 \,\mu\text{l}$ of the probe onto each pretreated specimen.
- Cover specimens with a 22 mm x 22 mm coverslip (avoid trapped bubbles) and seal the coverslip.

We recommend using rubber cement (e.g., Fixogum) for sealing.

- Place slides on a hot plate or hybridizer and denature specimens for 5 min at 79°C.
- Transfer slides to a humidity chamber and hybridize overnight at 37°C (e.g., in a hybridization oven).

It is essential that specimens do not dry out during the hybridization step.

Post-hybridization

Perform post-hybridization processing (washing, detection, counterstaining, mounting, microscopy) according to the instructions for use of the <u>Zyto Dot 2C CISH Implementation Kit</u>.

13. Interpretation of results

Using the <u>Zyto Dot 2C CISH Implementation Kit</u>, hybridization signals of Digoxigenin-labeled polynucleotides appear as dark green colored distinct dots (distal to the USP6 breakpoint region), and Dinitrophenyl-labeled polynucleotides appear as bright red colored distinct dots (proximal to the USP6 breakpoint region).

Normal situation: In interphases of normal cells or cells without a translocation involving the USP6 gene region, two red/green fusion signals appear (see Fig. 2).

Aberrant situation: One USP6 gene region affected by a translocation is indicated by one separate green signal and one separate red signal (see Fig. 2).

Overlapping signals may appear as brown signals.

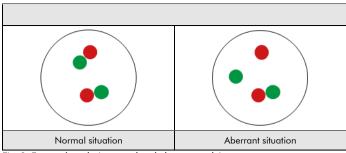


Fig. 2: Expected results in normal and aberrant nuclei

Genomic aberrations due to small deletions, duplications or inversions might result in inconspicuous signal patterns.

Other signal patterns than those described above may be observed in some abnormal samples. These unexpected signal patterns should be further investigated.

Please note:

- Due to decondensed chromatin, single CISH signals can appear as small signal clusters. Thus, two or three signals of the same size, separated by a distance ≤ 1 signal diameter, should be counted as one signal.
- Prior to signal enumeration, the specimen should be scanned for any possible intratumoral heterogeneity at 100- to 200-fold magnification.
- Visualization of signals should be performed at least at 400-fold magnification resulting in easily visible signals. A 630-fold magnification is recommended for probes detecting chromosomal breaks. Do not use contrast enhancing filter lenses as this might distort the signal color. To obtain signals in bright colors, open the aperture diaphragm. Be sure to focus up and down when evaluating a nucleus, as red and green signals might be located on top of each
- Do not evaluate areas of necrosis, overlapping nuclei, over-digested nuclei and nuclei with weak signal intensity.
- Due to mitosis, additional signals may be visible even in a small percentage of non-neoplastic cells. Occasionally, nuclei with missing signals may be observed in paraffin-embedded specimens due to cutting artefacts.
- A negative or unspecific result can be caused by multiple factors (see chapter 17. "Troubleshooting").
- In order to correctly interpret the results, the user must validate this product prior to use in diagnostic procedures according to national and/or international guidelines.

14. Recommended quality control procedures

In order to monitor correct performance of processed specimens and test reagents, each assay should be accompanied by internal and external controls. If internal and/or external controls fail to demonstrate appropriate staining, results with patient specimens must be considered invalid.

Internal control: Non-neoplastic cells within the specimen that exhibit normal signal pattern, e.g., fibroblasts.

External control: Validated positive and negative control specimens.

15. Performance characteristics

The performance of the probe was determined by comparison against the corresponding IVD approved FISH probe. The concordance was 100%.

Accuracy: The accuracy was calculated as 100%.

Analytical sensitivity: The analytical sensitivity was calculated as 100%.

Analytical specificity: The analytical specificity was calculated as 100%.

16. Disposal

The disposal of reagents must be carried out in accordance with local regulations.

17. Troubleshooting

Any deviation from the operating instructions can lead to inferior staining results or to no staining at all.

Weak signals or no signals at all		
Possible cause	Action	
Cell or tissue sample has not been properly fixed	Optimize fixing time and fixative	
Heat pretreatment, proteolysis, hybridization, denaturation, stringency wash or antibody-incubation temperature not correct	Check temperature of all technical devices used, using a calibrated thermometer. Use always the same number of slides in solutions with critical temperature	
Proteolytic pretreatment not carried out properly	Depending on multiple factors, e.g., nature and duration of fixing, thickness of sections, and nature of tissue/cells, different incubation times may be required. Ascertain the optimum time for pepsin incubation in pre-tests	
Hybridization time too short	Hybridize for at least 12 h; extend hybridization time if necessary	
Old dehydration solutions	Prepare fresh dehydration solutions	
Probe evaporation	When using a hybridizer, the use of the wet stripes/water filled tanks is mandatory. When using a hybridization oven, the use of a humidity chamber is required. In addition, the coverslip should be sealed completely, e.g., with Fixogum, to prevent drying-out of the sample during hybridization	
Incubation with chromogenic substrate too short	Extend incubation time	
Counterstaining time too long	The counterstaining time depends on the nature of the specimen and should be optimized accordingly. Avoid dark counterstaining, because it may obscure positive staining signals	
Bluing of counterstain not carried out properly	Use cold running tap water for bluing; do not use warm or hot water, or bluing reagents	

Signals too strong

Possible cause	Action
Proteolytic pretreatment carried out too long	Depending on multiple factors, e.g., nature and duration of fixing, thickness of sections, and nature of tissue/cells, different incubation times may be required. Ascertain the optimum time for pepsin incubation in pre-tests
AP-Red Solution incubation time not correct	If required, the incubation time can be shortened down to 5 min. Do not heat substrate solution over 25°C; incubate at room temperature only
HRP-Green solution incubation time not correct	If required, the incubation time can be shortened down to 7 min. Do not heat substrate solution over 25°C; incubate at room temperature only

Red signals too weak

rea signals lee weak		
Possible cause	Action	
AP-Red Solution was exposed to strong direct light	Prepare and use AP-Red Solution protected from strong direct light	
AP-Red Solution was prepared too early	Prepare prior to immediate use	
AP-Red Solution incubation time not correct	If required, the incubation time can be extended up to 15 min	
Insufficient preparation of chromogenic substrate	Do not increase volume of Solution A	

Green signals too weak

Green signals 100 weak		
Possible cause	Action	
Incubation time of any washing steps after staining with HRP-Green too long	Do not exceed given incubation times	
HRP-Green solution incubation time not correct	If required, the incubation time can be extended up to 15 min	
Insufficient preparation of chromogenic substrate	Do not increase volume of Solution A	

Signals fade or merge

Possible cause	Action
An unsuitable mounting solution has been used	Use only the mounting solution provided with the kit or xylene-based mounting solutions free of any impurities; do not use coverslip tape
Sections were not dehydrated properly	Use fresh ethanol and xylene solutions; use only xylene of "pure" quality

Uneven or in some parts only very light staining

Possible cause	Action
Incomplete dewaxing	Use fresh solutions; check duration of dewaxing times
Reagent volume too small	Ensure that the reagent volume is large enough to cover the tissue area
Air bubbles caught before hybridization or during mounting	Avoid air bubbles

Inconsistent results

Possible cause	Action
Insufficient drying before probe application	Extend air-drying
Too much water/wash buffer on tissue prior to application of pepsin, antibodies and/or color substrates	Ensure that excess liquid is removed from tissue section by blotting or shaking it off the slide. Small amounts of residual water/wash buffer do not interfere with the test
Variations in tissue fixation and embedding methods	Optimize fixation and embedding methods
Variations in tissue section thickness	Optimize sectioning

Marphalagy degraded

Morphology degraded		
Possible cause		Action
Cell or tissue sample l been properly fixed	has not	Optimize fixing time and fixative
Proteolytic pretreatme carried out properly	nt not	Optimize pepsin incubation time; decrease if necessary

Cross hybridization signals; noisy background

Possible cause	Action
Stringency wash temperature not correct	Check temperature of the technical devices used, using a calibrated thermometer. Use always the same number of slides in the jar. We recommend not to use more than eight slides per jar for heat incubation steps
Slides not thoroughly rinsed	Use fresh and sufficient wash buffer and deionized or distilled water where indicated
Sections dried out any time during or after hybridization	Avoid sections being dried out; use humidity chamber; seal coverslip properly

Prolonged substrate incubation time	Shorten substrate incubation time
Incomplete dewaxing	Use fresh solutions; check duration of dewaxing
Proteolytic pretreatment too strong	Optimize pepsin incubation time
Slides cooled to room temperature before hybridization	Transfer the slides quickly to hybridization temperature

Overlapping signals

Possible cause	Action
Inappropriate thickness of tissue sections	Prepare 3-5 μ m microtome sections

Specimen floats off the slide

operation near on the char	
Possible cause	Action
Unsuitable slide coating	Use appropriate (positively charged) slides
Proteolytic pretreatment too strong	Shorten pepsin incubation time

18. Literature

- Erickson-Johnson MR, et al. (2011) Lab Invest 91: 1427-33.
- Nakamura T, et al. (1988) Oncogene Res 2: 357-70.
- Oliveira AM, et al. (2004) Cancer Res 64: 1920-3.
- Oliveira AM, et al. (2005) Oncogene 24: 3419-26
- Wilkinson DG: In Situ Hybridization, A Practical Approach, Oxford University Press (1992) ISBN 0 19 963327 4.

19. Revision

Please refer to www.zytovision.com for the most recent instructions for use as well as for instructions for use in different languages.

Our experts are available to answer your questions. Please contact <u>helptech@zytovision.com</u>



ZytoVision GmbH Fischkai 1

27572 Bremerhaven/ Germany Phone: +49 471 4832-300 Fax: +49 471 4832-509 www.zytovision.com

Email: info@zytovision.com

Trademarks:

ZytoVision® and ZytoDot® are trademarks of ZytoVision GmbH.